


zepter[®]
INTERNATIONAL
LIVE BETTER • LIVE LONGER



PRO SKIN KIT

VB-PSK

by vacSy[®]

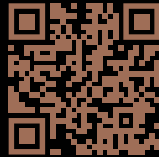

Congratulations on purchasing the ProSkin Kit!

Embark on a journey of self-care and indulge in an innovative skincare experience with the ProSkin Kit, a collection of beauty accessories meticulously designed to pamper and revitalize your skin.

The gentle suction created by vacuum therapy stimulates cellular regeneration, enhances blood circulation, and promotes collagen production, while the touch of real diamond powder brightens and rejuvenates.

With the ProSkin Kit, you achieve results that rival those of the most prestigious beauty salons, right in the comfort of your own home.

**PRO
SKIN
KIT** *by* vacSy[™]

CONTENTS

1. Introduction

2. Safety Instructions and Warnings

3. ProSkin Kit

- 3.1. Exfoliating Head
- 3.2. Blackheads Head
- 3.3. Face-lift Head
- 3.4. Cupping Head

4. Usage

5. Maintenance

6. Compliance declaration

The mechanism behind ProSkin

ProSkin employs a gentle suction mechanism that produces five healthy effects, ensuring a spa-like experience in the comfort of your home.

Each of the accessories from the SkinPro Kit is meticulously crafted to cater to specific beauty and wellness needs, ensuring a comprehensive and personalized experience.



For All Skin
Types



Renews
Skin Cells



Smoother
and Velvety
Skin

**PRO
SKIN
KIT** *by* vacSy[®]

About this instruction manual

Please read the safety instructions carefully and only use the product as described in these instructions to avoid accidental injury or damage.



This symbol warns you of the risk of injury.



This symbol warns you of the risk of injury from electricity.

The package contains:

- exfoliating head
- blackheads head
- facelift head
- 2 cupping heads
- 25 replaceable filters
- 1 accessories pouch

1. Introduction

Welcome to ProSkin Kit - Where beauty meets technology.








The first-ever diamond powder microdermabrasion device cleanses pores, eliminates blackheads, combats the signs of aging, stimulates collagen and elastin production, promotes cell renewal, and enhances skin radiance. A mere 10 minutes of weekly facial massage will transform your skin in just 4 weeks, leaving it healthier and more radiant.





It is important to mention that these ProSkin Kit accessories can only be used in combination with the Vacsy Vacuum Pump Deluxe, VB-P from Zepter.

The pump is equipped with a sensor that will automatically recognize the accessories and activate the various treatment functions.

2. Safety Instructions and Warnings

-  Read the manual before use and keep it for future reference.
-  The device is not a toy, keep out of reach of children.
-  Do not use if the skin is swollen, irritated or inflamed.
-  Make sure that the polymer filter is positioned correctly, the filter, hands, and skin must be dry. The skin may become slightly red after use. If redness persists, consult a doctor about further usage.
-  Do not use peeling, abrasive or other irritating products before and after using the ProSkin Kit.
-  Not suitable for pregnant women; children; people suffering from epilepsy, tumors, or inflammation; people with a pacemaker; people sensitive to heat and vibration; persons suffering from eczema or dermatological diseases.
-  Limit usage to a maximum of twice a week.

During usage

-  **Avoid abrupt lifting from the skin.**
-  **Do not stop in the same area for more than 2 seconds.**

3. ProSkin Kit

3.1. Exfoliating Head (Diamond Dust Surface): Reveal your radiance



Aspiration power:



Toxins removal:



Microcirculation improvement:



It features a microcrystalline surface made of real diamond dust, offering numerous benefits by gently removing dull, dead skin cells, setting the perfect foundation for your skincare routine.

Recommended Usage

Suitable for all skin types and can be used all over the face. The recommended duration for each use is 8-10 minutes for optimal results.



! Move the device back and forth on the skin.

! Use slow and fluid movements.



**PRO
SKIN
KIT** *by* **vacSy**

3.2. Blackheads Head: Banish blemishes, welcome clarity



Aspiration power:
████████████████████

Toxins removal:
██████████████████

Microcirculation improvement:
██████████████████

Designed to remove blackheads, reduce sebum production, and deeply clean pores. Embrace the confidence that comes with clear, radiant, and beautiful skin.

Recommended Usage

Ideal for oily and mixed skin types, this head is specifically recommended for use on the nose area, due to its potency.



3.3. Face-lift Head: Sculpt and rejuvenate



Aspiration power:
 ██████████

Toxins removal:
 ██████████

Microcirculation improvement:
 ██████████

Its unique shape is perfect for massaging cutaneous tissues and activating blood and lymphatic microcirculation. It aids in collagen and elastin production, reducing wrinkles and having a lifting effect.

Recommended Usage

Suitable for use on cheeks, cheekbones, forehead, and jawline.



! Move the device back and forth on the skin.

! Use slow and fluid movements.



**PRO
SKIN
KIT** *by* **vacSy**

3.4. Cupping Head: Boost your overall well-being



Aspiration power:
████████████████████

Toxins removal:
██████████████████

Microcirculation improvement:
████████████████████


Introduction to Cupping

Cupping is an ancient oriental technique using cups to create a void, promoting collagen and elastin production. It's effective for beauty treatments, reducing cellulite, and aiding athletes post-training by:

- Promoting elastin and collagen production
- Reducing oedemas, limb, and muscle pain
- Treating cellulite and stretch marks


Recommended Usage

The Cupping Head can be used on arms, shoulders, back, abdomen, hips, thighs, and glutes.

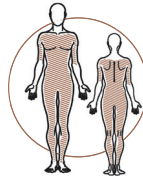

Cupping Treatments

Before starting your treatment, remember to nourish your skin with an oil or a body cream. This will ensure that the cupping head glides effortlessly across your skin, minimizing any discomfort or irritation.

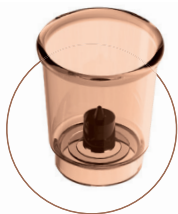

Traditional Mode

- ① Place the cup on the affected area.
- ① Create a void, then remove the vacuum.
- ① Leave the cup on the skin until the void fades.
- ① This mode is perfect to treat muscle and limb pain.

Massage Mode

- ① Follow the instructions for Traditional Mode.
- ① Manually move the head for approx. 10 minutes.
- ① Ideal for cellulite, stretch marks, and relaxing massage.



Remedy for
Pain



Reducing
Cellulite




Body
Cupping




Stimulates
Circulation

Caution and Suitability

Avoid using on bony or very thin individuals and areas with wounds or varicose veins.

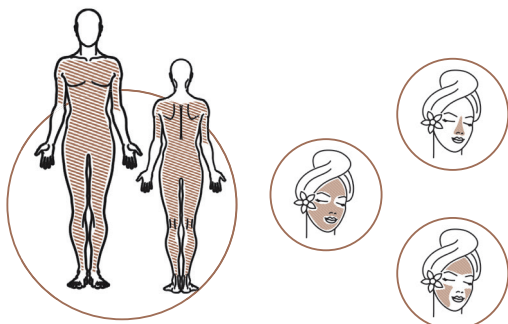

4. Usage

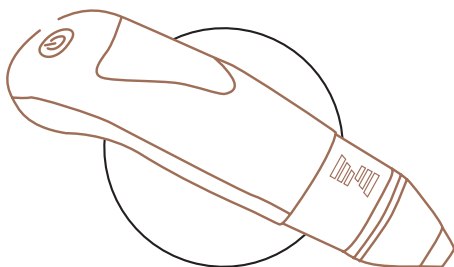

a) Pre-treatment Recommendations

For optimal results and safety:

- ① Fully charge the Vacuum pump before the first use.
 - ① Use a low suction level for initial treatments.
 - ① Open pores with steam or a warm towel before each treatment.
-

b) Operation Steps

- ① Choose the head suitable for your treatment.
 - ① Insert the filter inside the head.
 - ① Connect the head to the Vacsy Pump.
 - ① Click one time on the power button to set the suction speed at level 1.
 - ① Each subsequent click increases the speed up to level 4.
 - ① Hold down the power button for 5 seconds to turn off the pump.
-

**PRO
SKIN
KIT** *by* **vacSy**

c) Summarized Steps for Perfect Skin

- ① Use makeup remover and wash your face.
- ① Open pores with steam.
- ① Begin ProSkin Kit treatment.
- ① Apply face cream or nourishing oil.
- ① Enjoy your radiant skin! Repeat up to twice a week.

d) Frequency and Results

- ① Results will appear in 4-6 weeks of consistent use.
Follow the guidelines to unleash your skin's natural glow.
-

5. Maintenance

ProSkin Kit Cleaning

- ① Wash and wipe heads after each use.
-

Vacuum Pump Cleaning

- ① Use a wet towel to clean the Vacuum Pump.
 - ① Do not soak the Vacuum Pump; it may cause irreparable damage.
-

Filters Replacement

- ① Change filters after 6-8 uses.

6. Declaration of compliance

Zepter Group hereby declares that all chemical substances utilized in the manufacturing of these accessories are in full compliance with all relevant rules and regulations set forth by the European Union.

This declaration attests to our commitment to meeting the highest standards of quality and safety in accordance with applicable EU directives.

Discover the Zepter Healthy Lifestyle products



BIOPTRON⁺
HYPERLIGHT THERAPY SYSTEM

MedAll

SAFE AND EFFECTIVE
MEDICAL TREATMENT
FOR HEALTH,
BEAUTY & WELL-BEING




**SWISSO
LOGICAL**

SO SWISS, SO LOGICAL.



**LA
DANZA**

HYPER-RETIJUVENATE

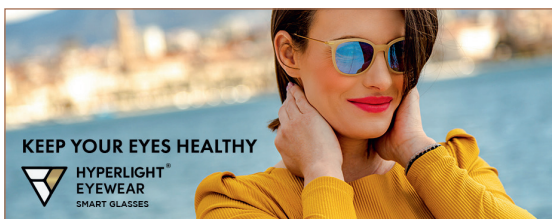

VB-PSK

PRO SKIN KIT

by vacSy®

zepter®
INTERNATIONAL
LIVE BETTER • LIVE LONGER
www.zepter.com

Home Art & Sales Services AG
Sihleggstrasse 23
8832 Wollerau
Switzerland